



**Mess.-Nr.** 60284 **Datum** Dienstag, 29 Januar 2013 04:49 **Bearbeitungs- Nr.:**  
**Material** NaCl (trocken) **PI 5cl 400rpm 500ml ZrO2 20mm** **M130007**  
**Bemerkung** predis. in IPA 1 min Ultrasonic

**Berechnung** Fraunhofer broad **TradeOff** 1.000000e+003  
**Scans** 100 **Scans** 100 **Messkanäle** 108 **Strahlabsorption** 17.00 %  
**Feinbereich** **Grobbereich**  
**Messbereich** 0,01 - 2000.00 [ $\mu\text{m}$ ] **Pumpe** 5 **Ultraschall** 0

| %  | < $\mu\text{m}$ |
|----|-----------------|
| 5  | 0.88            |
| 10 | 1.24            |
| 15 | 1.61            |
| 20 | 2.04            |
| 25 | 2.56            |
| 30 | 3.24            |
| 40 | 5.38            |
| 50 | 9.11            |
| 60 | 14.37           |
| 70 | 20.50           |
| 80 | 27.56           |
| 90 | 36.91           |
| 95 | 44.48           |
| 98 | 52.87           |
| 99 | 58.34           |

| %     | < $\mu\text{m}$ |
|-------|-----------------|
| 0.1   | 0.1             |
| 1.0   | 0.5             |
| 6.7   | 1.0             |
| 38.6  | 5.0             |
| 51.9  | 10.0            |
| 61.1  | 15.0            |
| 69.2  | 20.0            |
| 76.6  | 25.0            |
| 83.0  | 30.0            |
| 88.3  | 35.0            |
| 92.4  | 40.0            |
| 97.2  | 50.0            |
| 99.2  | 60.0            |
| 99.9  | 70.0            |
| 100.0 | 80.0            |

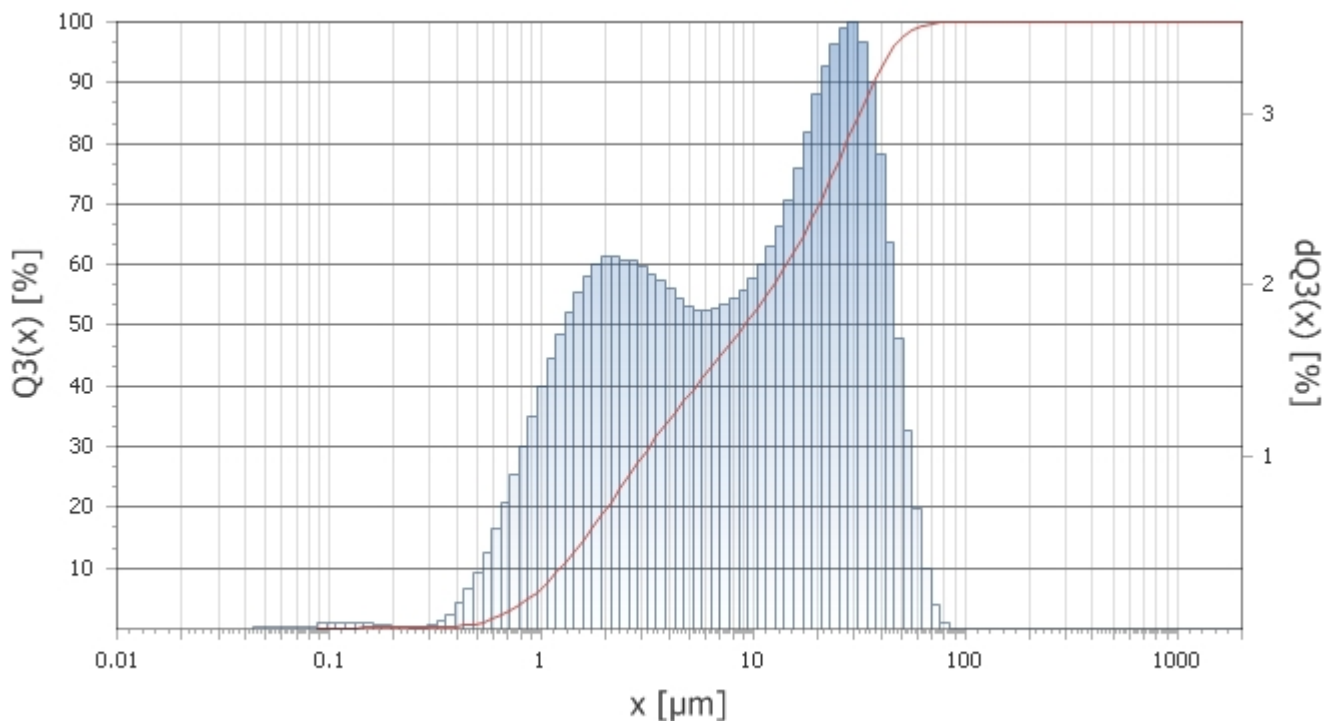

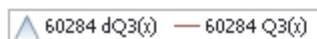
 60284 dQ3(x) — 60284 Q3(x)