


Mess.-Nr. 60293 **Datum** Mittwoch, 30 Januar 2013 08:50 **Bearbeitungs- Nr.:**
Material NaCl (trocken) **PI 5cl 400rpm 500ml ZrO2 20mm** **M130007**
Bemerkung predis. in IPA 1 min Ultrasonic

Berechnung Fraunhofer broad **TradeOff** 1.000000e+003
Scans 100 **Scans** 100 **Messkanäle** 108 **Strahlabsorption** 16.00 %
Feinbereich **Grobereich**
Messbereich 0,01 - 2000.00 [μm] **Pumpe** 5 **Ultraschall** 0

%	< μm
5	0.75
10	1.02
15	1.27
20	1.54
25	1.84
30	2.18
40	3.10
50	4.68
60	7.57
70	11.64
80	16.85
90	24.65
95	31.37
98	38.99
99	43.97

%	< μm
0.2	0.1
1.6	0.5
9.6	1.0
27.5	2.0
51.4	5.0
66.2	10.0
76.8	15.0
84.7	20.0
90.3	25.0
94.2	30.0
96.7	35.0
98.2	40.0
99.2	45.0
99.6	50.0
100.0	60.0

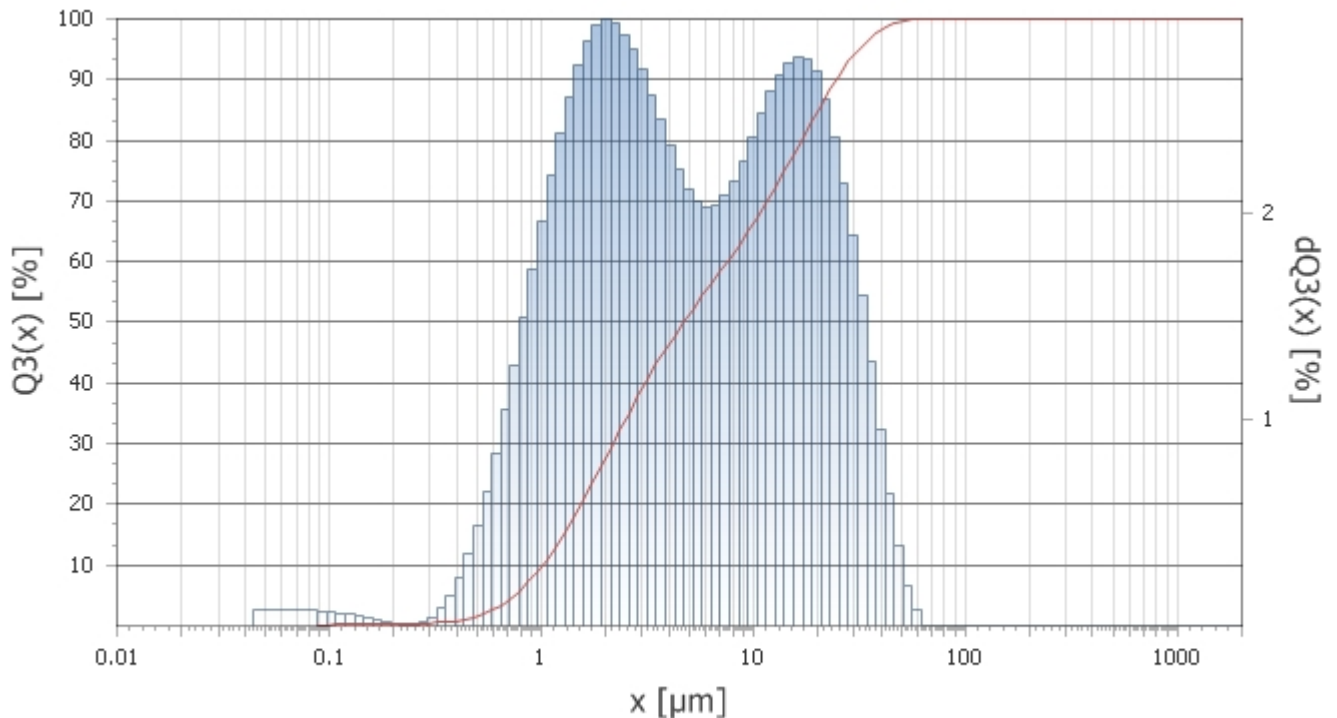

▲ 60293 dQ3(x) — 60293 Q3(x)