


Mess.-Nr.	60300	Datum	Mittwoch, 30 Januar 2013 10:38	Bearbeitungs- Nr.:
Material	NaCl		P5cl 500ml 20mm ZrO2 30min (10' SBÖ)	M130007
Bemerkung	predisp. in IPA 1 min Ultrasonic			

Berechnung	Fraunhofer	broad	TradeOff	1.000000e+003
Scans				
Feinbereich	100	Scans Grobbereich	100	Messkanäle 108
Messbereich	0,01	- 2000.00 [µm]	Pumpe	5
			Ultraschall	0

%	< µm
5	0.72
10	0.97
15	1.20
20	1.43
25	1.69
30	1.97
40	2.69
50	3.78
60	5.73
70	9.37
80	14.90
90	22.78
95	29.23
98	36.21
99	40.60

%	< µm
0.2	0.1
1.7	0.5
10.6	1.0
30.5	2.0
57.0	5.0
71.3	10.0
80.2	15.0
87.0	20.0
92.0	25.0
95.4	30.0
97.6	35.0
98.9	40.0
99.6	45.0
99.9	50.0
100.0	60.0

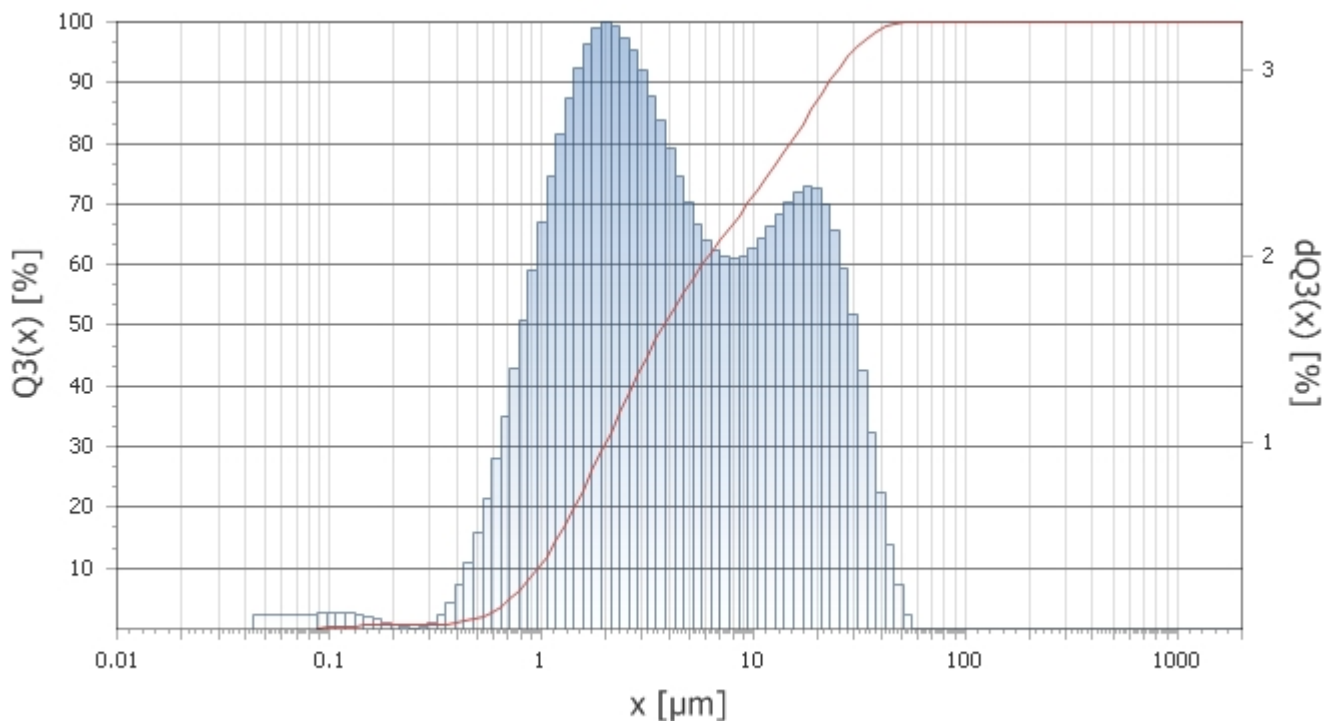

▲ 60300 dQ3(x) — 60300 Q3(x)