


Mess.-Nr. 60272 **Datum** Dienstag, 29 Januar 2013 03:50 **Bearbeitungs- Nr.:**
Material NaCl (trocken) **PI 5cl 400rpm 500ml ZrO2 30mm** **M130007**
Bemerkung predis. in IPA 1 min Ultrasonic

Berechnung Fraunhofer broad **TradeOff** 1.000000e+003
Scans 100 **Scans** 100 **Messkanäle** 108 **Strahlabsorption** 16.00 %
Feinbereich **Grobbereich**
Messbereich 0,01 - 2000.00 [µm] **Pumpe** 5 **Ultraschall** 0

%	< µm
5	0.90
10	1.46
15	2.20
20	3.12
25	4.10
30	5.08
40	7.14
50	9.54
60	12.69
70	17.29
80	24.08
90	34.19
95	42.37
98	51.21
99	56.95

%	< µm
0.0	0.1
1.2	0.5
5.9	1.0
29.6	5.0
51.6	10.0
65.5	15.0
74.5	20.0
81.1	25.0
86.4	30.0
90.6	35.0
93.8	40.0
97.7	50.0
99.3	60.0
99.9	70.0
100.0	80.0

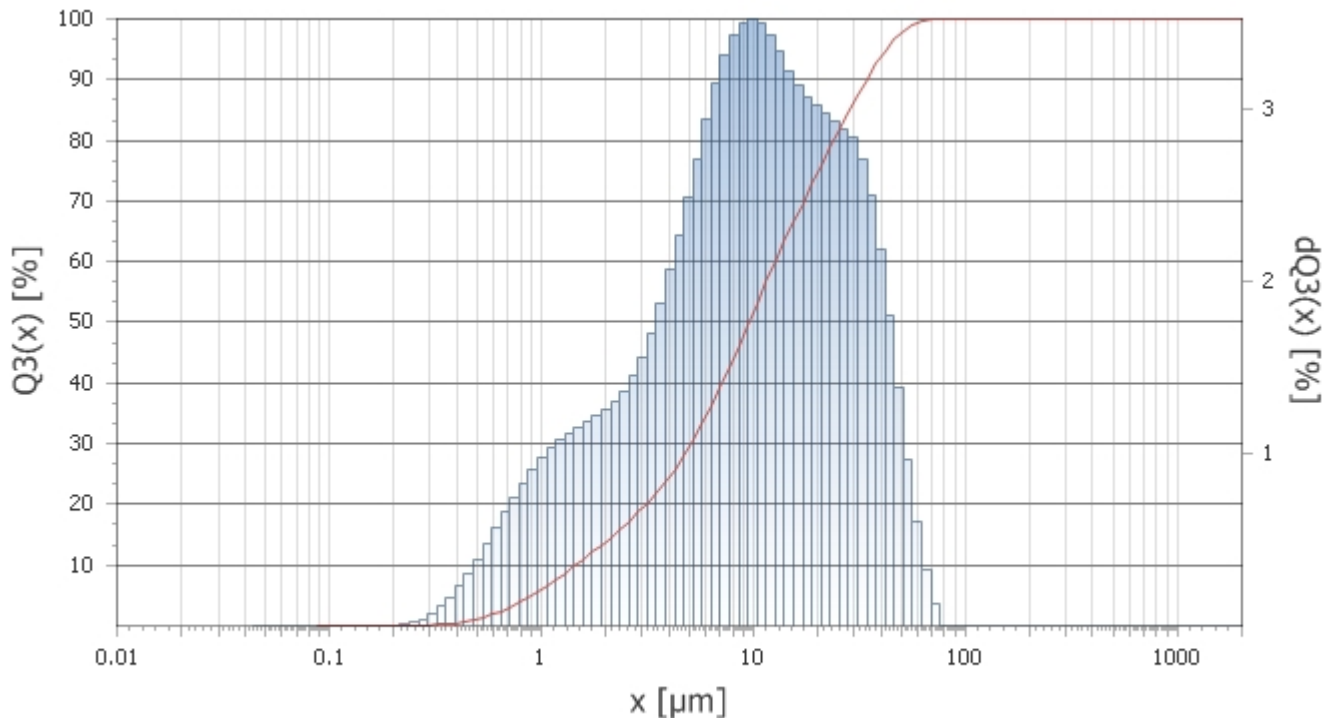

▲ 60272 dQ3(x) — 60272 Q3(x)