


Mess.-Nr. 60275 **Datum** Dienstag, 29 Januar 2013 03:55 **Bearbeitungs- Nr.:**
Material NaCl (trocken) **PI 5cl 400rpm 500ml ZrO2 30mm** **M130007**
Bemerkung predis. in IPA 1 min Ultrasonic

Berechnung Fraunhofer broad **TradeOff** 1.000000e+003
Scans 100 **Scans** 100 **Messkanäle** 108 **Strahlabsorption** 18.00 %
Feinbereich **Grobbereich**
Messbereich 0,01 - 2000.00 [μm] **Pumpe** 5 **Ultraschall** 0

%	< μm
5	0.80
10	1.15
15	1.46
20	1.79
25	2.14
30	2.55
40	3.62
50	5.48
60	9.24
70	14.50
80	20.48
90	28.81
95	35.77
98	43.54
99	48.72

%	< μm
0.5	0.1
2.0	0.5
7.8	1.0
48.0	5.0
61.5	10.0
70.9	15.0
79.3	20.0
86.0	25.0
91.0	30.0
94.6	35.0
96.9	40.0
99.2	50.0
99.8	60.0
100.0	70.0
100.0	80.0

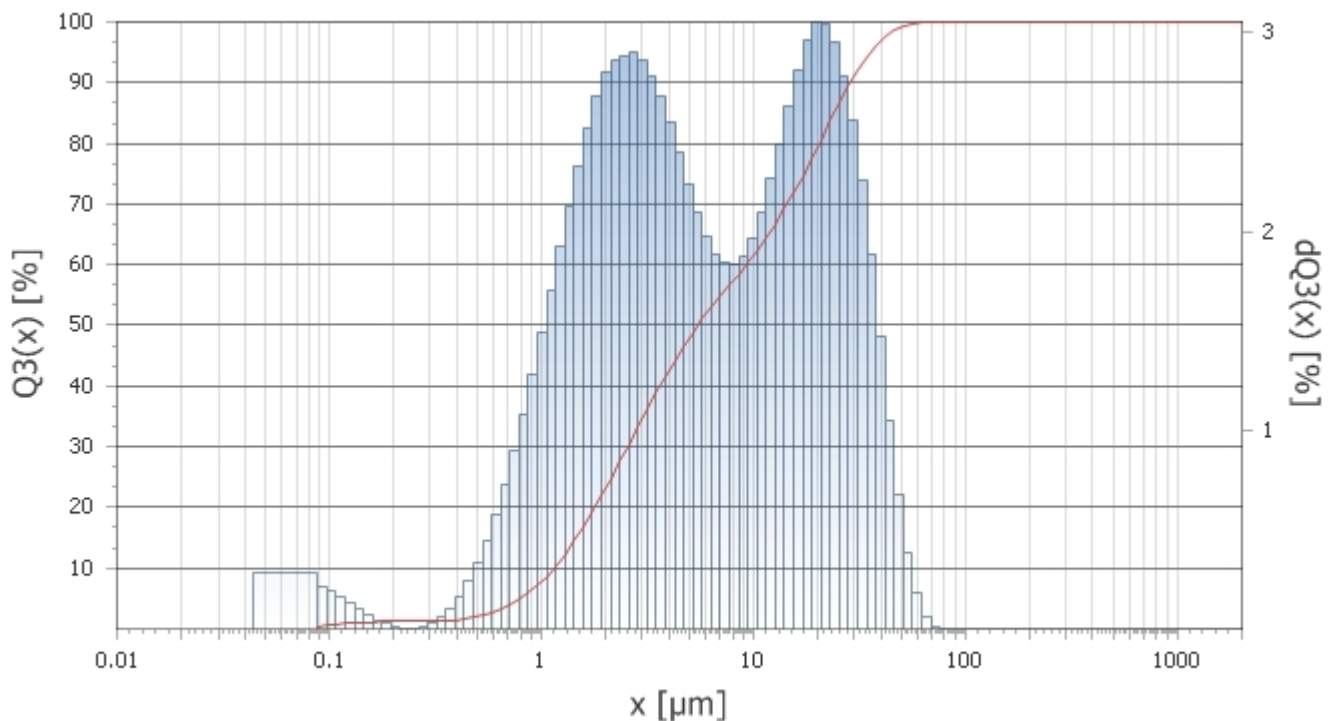

 60275 dQ3(x) — 60275 Q3(x)