


**Meas.Nr.** 62689      **Date** 04.03.2013 13:14      **Personal Report No.:**  
**Material** Zircon B1      60min P7pl 0,1mm ZrO2      **M130032**  
**Description** vordisp. in 0,1% Na4P2O7 mit 1min Ultraschall LABORETTE 17

**Calculation** Mie      n = 1.98970      -i 0      **broad TradeOff** 1.000000e+002  
**Scans Fine** 100      **Scans Coarse** 0      **Channels** 108      **Beam Obscuration** 15.00 %  
**Meas. Range** 0.01 [µm] - 42.30 [µm]      **Pump** 4      **Ultrasonics** 10

%	< µm
5	0.144
10	0.161
15	0.176
20	0.190
25	0.204
30	0.220
40	0.258
50	0.335
60	0.445
70	0.538
80	0.636
90	0.764
95	0.875
98	1.002
99	1.082

%	< µm
0.0	0.08
0.1	0.10
6.5	0.15
23.6	0.20
46.3	0.30
55.7	0.40
65.8	0.50
76.5	0.60
85.5	0.70
91.8	0.80
95.7	0.90
98.0	1.00
99.7	1.20
100.0	1.50
100.0	2.00

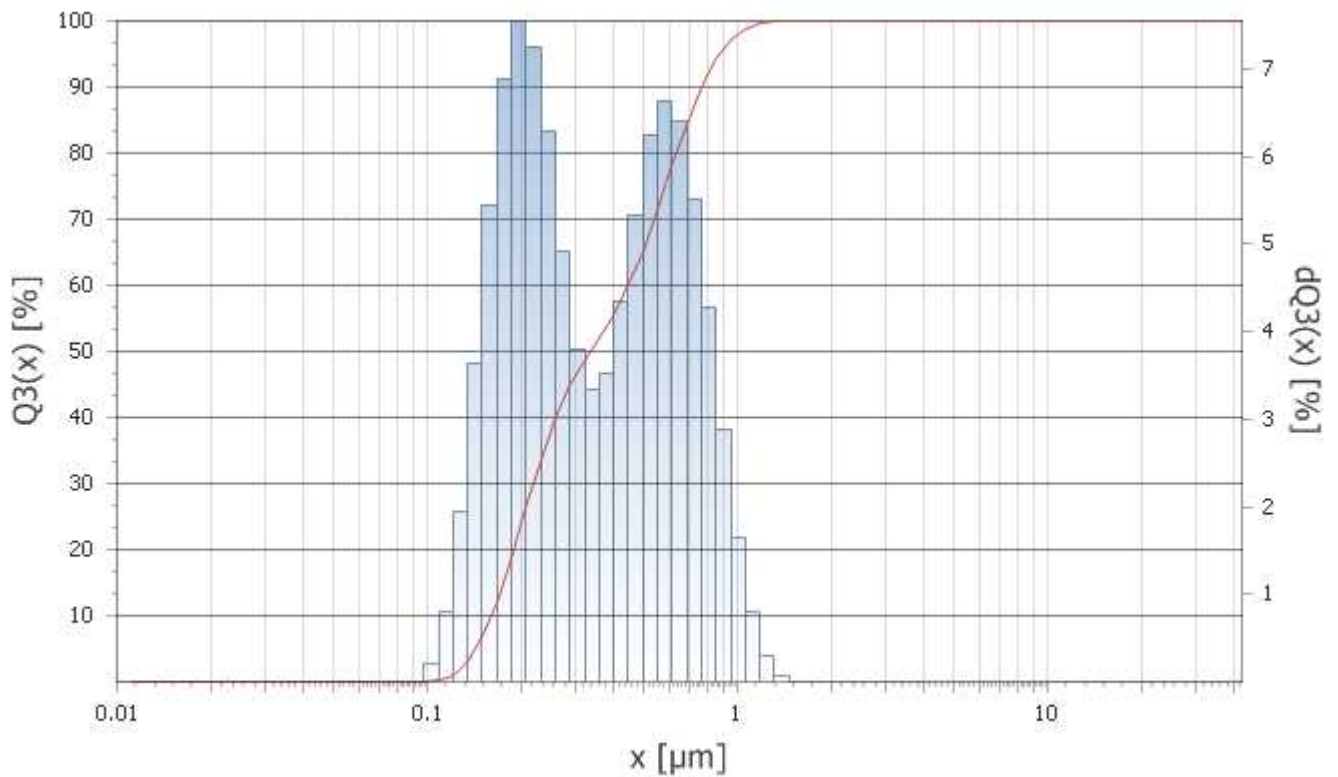

▲ 62689 dQ3(x) — 62689 Q3(x)